


Sampling Instructions – ALPACA

generatio

Animal Trust Center
Heidelberg - Germany
www.generatio.com

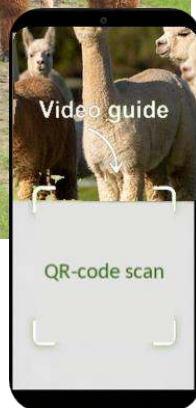

- ✓ **Hair roots**
- ✓ **EDTA blood sample (vet. staff only)**



Sampling Instructions - ALPACA

generatio

Animal Trust Center
Heidelberg - Germany
www.generatio.com

Hair roots

Important notes:

1. For officially recognised tests, sampling must be performed by a veterinarian.
2. When submitting samples for a breeding organisation, follow the applicable regulations.


Preparation of the Hair Sample Envelope

Use paper envelope of your choice

Write required information outside of the envelope

Required information:

- ☐ Animal name
- ☐ Phone for inquiries
- ☐ Order number of online order

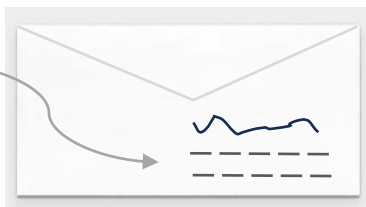

Generatio GmbH
Blumenstraße 49
69115 Heidelberg
Germany

Contact
support@generatio.com
+49 (0) 6221 389 353 0 Tel
+49 (0) 6221 389 353 1 Fax

Sampling Instructions - ALPACA

Hair roots

Important notes:

1. For officially recognised tests, sampling must be performed by a veterinarian.
2. When submitting samples for a breeding organisation, follow the applicable regulations.

Sampling

➡ Obtain two finger-thick hair bundles with roots from the tail area

➡ Hold hair bundles together with adhesive strip

Video
instructions



Notes

👉 Requirements for hair before sampling:

- Dry
- Free of styling products

Completed sample envelope with camel hairs

Paper envelope labelled on the outside:

- ✓ Animal name
- ✓ Phone for inquiries
- ✓ Order number of online order

In the envelope:

- ✓ two 'finger-thick' bundles
- ✓ hairs with hair roots bundled with adhesive strip



If you are submitting hair samples from multiple animals, please place each sample in a separate envelope!

Sampling Instructions - ALPACA

generatio


Animal Trust Center
Heidelberg - Germany
www.generatio.com

Hair roots

Important notes:

1. For officially recognised tests, sampling must be performed by a veterinarian.
2. When submitting samples for a breeding organisation, follow the applicable regulations.

Prepare sample shipment

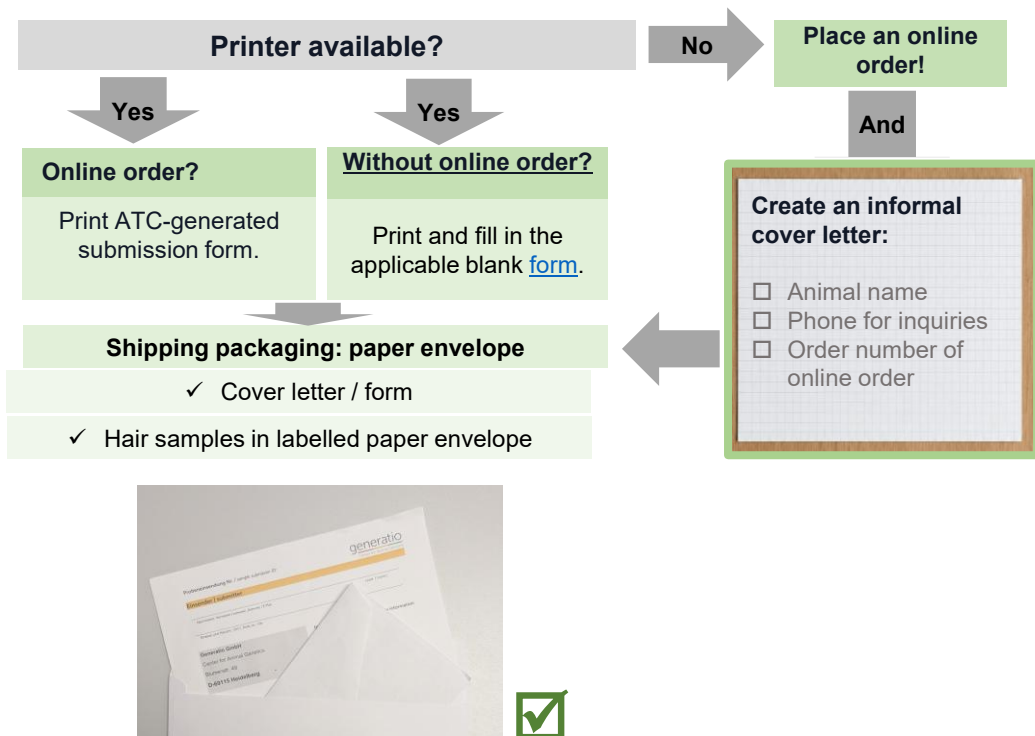

Shipping instructions



Generatio GmbH
Blumenstraße 49
69115 Heidelberg
Germany

Contact
support@generatio.com
+49 (0) 6221 389 353 0 Tel
+49 (0) 6221 389 353 1 Fax

Sampling Instructions - ALPACA


EDTA blood sample (vet. personnel only)

Important notes:

1. *For officially recognised tests, sampling must be performed by a veterinarian.*
2. *When submitting samples for a breeding organisation, follow the applicable regulations.*

Preparation

☐ Label tube for EDTA blood

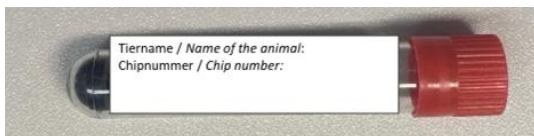

Sampling

- ➡ Draw 2–3 ml of blood.
- ➡ Carefully mix blood with EDTA by inverting several times.
- ➡ Please do not additionally seal the EDTA blood tube with Parafilm.
- 👉 Protect sample from direct sunlight.
- 👉 If shipping later than 2 days: store sample tube in the refrigerator.

Sampling Instructions - ALPACA

generatio

Animal Trust Center
Heidelberg - Germany
www.generatio.com

EDTA blood sample (vet. personnel only)

Important notes:

1. For officially recognised tests, sampling must be performed by a veterinarian.
2. When submitting samples for a breeding organisation, follow the applicable regulations.

Completed EDTA blood tube



Pack the labelled EDTA blood in a break-resistant transport sleeve.

Prepare sample shipment

Printer available?

No

Place an online order!

Yes

Yes

Online order?

Print ATC-generated submission form.

Without an online order?


Print and fill in the applicable blank [form](#).

And

Create an informal cover letter:

- ☐ Animal name
- ☐ Phone for inquiries
- ☐ Order number of online order

Packaging: tear-resistant / padded envelope

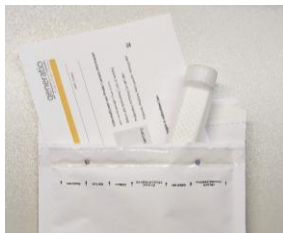

Shipping instructions



Generatio GmbH
Blumenstraße 49
69115 Heidelberg
Germany

Contact
support@generatio.com
+49 (0) 6221 389 353 0 Tel
+49 (0) 6221 389 353 1 Fax